
## VITALOGRAPH Eco BVF BACTERIAL VIRAL FILTERS

Safer Pulmonary Function Testing During COVID-19 and Beyond



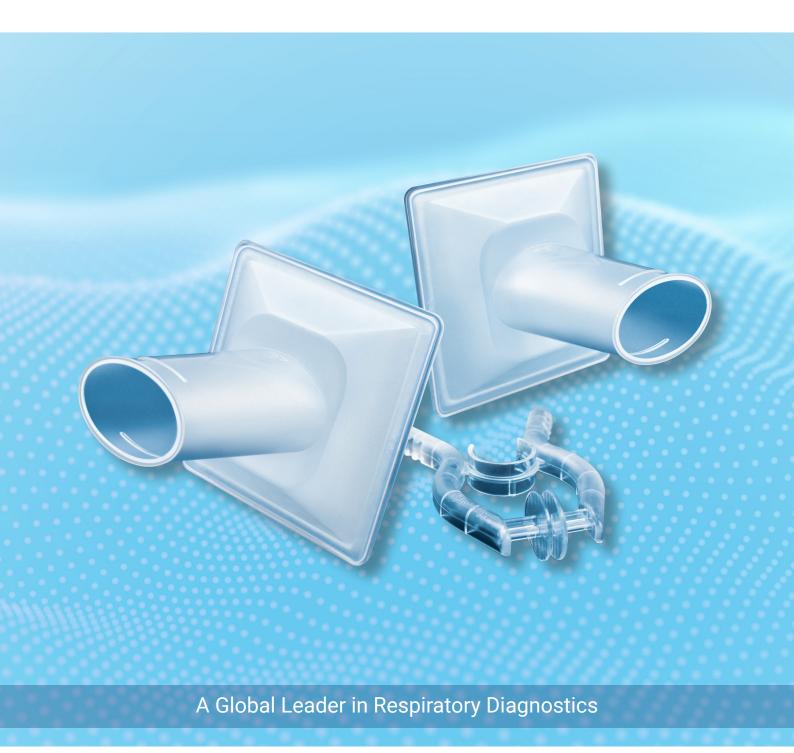


Vitalograph Eco BVF Bacterial Viral Filters

**High Performance Cross Contamination Protection** for Pulmonary Function Testing (PFT)

The Vitalograph® range of Eco BVFs combines independently validated >99.999% bacterial and viral cross-contamination efficiency for medical personnel, patients and equipment with a reduced environmental impact.

The range includes a 'Bite-Lip' design offering a more convenient alternative to conventional Filter/Bite-On mouthpiece combinations.

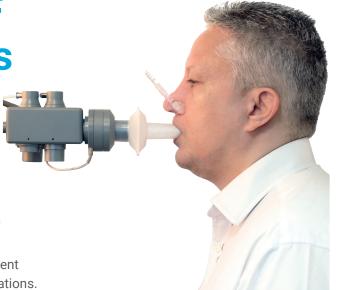


### **Reduced Environmental Impact**

- Single-material, 100% recyclable\* 🖧
- 80% reduction in energy expended during manufacturing\*\*
- 54% reduction in filter waste volume\*\*

# **Eco BVF Bite-Lip Option Reduces Consumable Costs**

Extended patient side with plastic bite-on lip reduces costs by eliminating the need for additional silicone

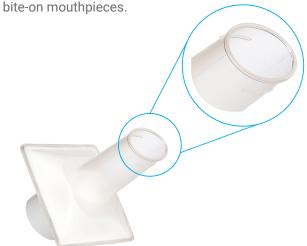

- \* Eco-BVF outer casing is Polypropylene and may be recycled with #5 materials. Eco-BVF filter is non-woven Polypropylene with acrylic strands. Check your local recycling capability to ensure this form of Polypropylene is accepted for recycling with #5 materials.
- \*\* compared to standard Vitalograph BVFs
- \*\*\* Thermoplastic Elastomer

### **Benefits**

- · Optimal protection for patients, staff and equipment
- Bite-Lip and standard Eco BVFs available for all leading makes of PFT equipment
- Filter dead space 55ml (excluding tapers)
- · Low flow impedance < 0.04kPa at 1L/s

## **Full Range of PFT Kits Available**

- PFT Kits available for Eco Bite-Lip BVFs and standard Eco BVFs
- Kits include choice of Eco BVF™ Bite-Lip + Nose Clips, or standard Eco BVF + TPE\*\*\* Bite On Mouthpiece + Nose Clips




Kits available for Jaeger/Carefusion/Vyaire, MediSoft, Ganshorn etc.



## **Cross-Contamination Efficiency for Bacteria and Viruses**

Vitalograph Eco BVFs filter droplets, aerosols, bacteria and viruses (including MRSA, TB, Influenza & SARS-CoV-2) out of the airstream with >99.999% cross-contamination efficiency.

Cardboard mouthpieces, disposable turbines and ultrasonic flow sensors do not prevent potentially germ carrying droplets and particulates from being expelled into the testing environment.

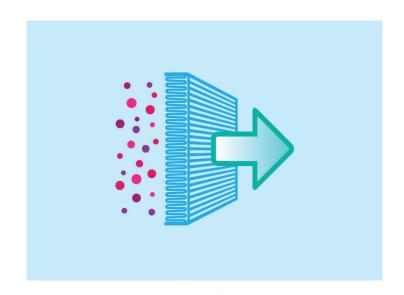

# Measurement Accuracy Assured for all Spirometer & PFT System Brands

You can be confident that the use of Eco BVFs will not give rise to measurement errors.

- Resistance (Flow Impedance) of the filter: Vitalograph filters are well below the maximum resistance value recommended by international ATS/ERS/ISO Spirometry Standards.
- 2. The Vitalograph Eco BVF has a comparatively small effective deadspace of 55ml to improve the validity and consistency of testing and minimise rebreathing.

## **Compliance with International Standards**

Vitalograph is a world leading manufacture of equipment for pulmonary function professionals. Vitalograph filters are CE marked and are designed and manufactured in line with the following requirements:

ISO 13485, FDA 21 CFR 820, CMDR, JPAL

Eco BVF and Eco BVF/Nose Clip Kits available for most office spirometers including; MIR, MicroMedical, Cosmed etc.



## **Technical Specifications**

Materials:PolypropylenePolypropyleneDimensions:68mm x 72mm68mm x 95mm

**Eco BVF** 

Weight: 11gm 25gm

Dead Space:55ml (exc. tapers)55ml (exc. tapers)Resistance:0.04kPa/(L/s) at 1L/s0.04kPa/(L/s) at 1L/s

**Cross-contamination efficiency:** >99.999% >99.999%

## **Ordering Information**

Eco BVF for PFT Equipment and Office Spirometers								
Product Model Number	Customer Order Number	Product Description	Product Box Size	Product Fits				
2820	28557	Eco BVF with Plastic Bite Lip	75	Carefusion / Jaeger / Vyaire				
	28550	Eco BVF with Plastic Bite Lip + Disposable Nose Clip	75	Carefusion / Jaeger / Vyaire				
	28552	Eco BVF + TPE* Bite-On Mouthpiece + Disposable Nose Clip	60	Carefusion / Jaeger / Vyaire				
	28554	Eco BVF with Plastic Bite Lip	75	Vitalograph / MediSoft / Office Spirometers				
	28553	Eco BVF with Plastic Bite Lip + Disposable Nose Clip	75	Vitalograph / MediSoft / Office Spirometers				
	28551	Eco BVF + TPE* Bite-On Mouthpiece + Disposable Nose Clip	60	Vitalograph / MediSoft / Office Spirometers				

**Eco BVF Bite-Lip** 

Eco BVF for Office Spirometers							
Product Model Number	Customer Order Number	Product Description	Product Box Size	Product Fits			
2820	28501	Eco BVF	100	Vitalograph / Cosmed / MicroMedical / MIR			
	28572	Eco BVF + Disposable Nose Clip	80	Vitalograph / MediSoft / Office Spirometers			



Adaptors							
Product Model Number	Customer Order Number	Product Description	Product Box Size	Product Fits			
2820	29748	Adaptor for Eco BVF, 28501, 28551, 28553, 28554, 28572	1	Koko			
	29808	Adaptor for Eco BVF, 28501, 28551, 28553, 28554, 28572	1	nSpire / ZAN			
	29828	Adaptor for Eco BVF, 28550, 28552, 28557	1	nSpire / ZAN			

<sup>\*</sup>Thermoplastic Elastomer

All other trademarks are the property of their respective owners.

56018\_5

## our services cover the globe

#### Vitalograph Inc.

Toll Free: 800 255 6626
Tel: (913) 730 3200
Fax: (913) 730 3232
contact@vitalograph.com
13310 West 99th Street
Lexexa, Kansas, 66215,
IISA

#### Vitalograph Ltd, Hong Kong/China

sales@vitalograph.cn P.O. Box 812 Shatin Central Post Office Hong Kong

#### Vitalograph (Ireland) Ltd

Tel: +353 65 686 4100 Fax: +353 65 682 9289 sales@vitalograph.ie Gort Road Business Park Ennis, Co Clare, V95 HFT4, Ireland

#### Vitalograph Ltd, International

Tel: +44 1280 827120 Fax: +44 1280 823302 sales@vitalograph.co.uk Maids Moreton, Buckingham MK18 1SW, England

#### Authorised Dealer Stamp

#### Vitalograph GmbH

Tel: +49 40 547391-0 Fax: +49 40 547391-40 info@vitalograph.de Rellinger Straße 64a D-20257 Hamburg Germany

### Vitalograph Ltd, UK

Tel: 01280 827110 Fax: 01280 823302 sales@vitalograph.co.uk Maids Moreton, Buckingham MK18 1SW, England

